

# StatSensor<sup>®</sup> CREAT

## Point-of-Care Whole Blood Creatinine and eGFR Testing



StatSensor Xpress



StatSensor Hospital Meter

Single-Use Creatinine Biosensor Test Strip

Virtually Painless Finger Stick Capillary Blood Sample, 1.2 µL

Results in 30 Seconds

Wide Measurement Range, 27-1056 µmol/L (0.3-12 mg/dL)

Calculated eGFR by CKD-EPI and MDRD Equations

*nova*<sup>®</sup>  
biomedical

# StatSensor® Whole Blood Creatinine and eGFR Testing



StatSensor® Creatinine Meter

## Nova StatSensor® Creatinine Measuring System

In response to the growing need to quickly and easily assess kidney function in many point-of-care (POC) settings, Nova Biomedical introduces StatSensor Creatinine, a handheld analyser and miniaturised, disposable biosensor for whole blood creatinine testing. StatSensor Creatinine incorporates patented Multi-Well™ technology adapted from Nova's hospital glucose monitoring system. StatSensor advanced technology enables simple, rapid, and accurate assessment of renal function by finger stick capillary blood sampling at the point of care.

### Features:

- Correlates with laboratory enzymatic creatinine methods
- Estimated glomerular filtration rate (eGFR) by CKD-EPI, MDRD and Cockcroft-Gault equations
- Simple test procedure (single-use, precalibrated disposables)
- Virtually painless, finger stick capillary blood sampling (1.2 µL)
- Results in 30 seconds
- Wide measurement range: 27 µmol/L to 1056 µmol/L (0.3-12 mg/dL)

## Creatinine with eGFR Enhances Assessment of Renal Function

Whole blood creatinine concentration has limitations when used alone to assess kidney function. While creatinine is the best single blood test for kidney function, it is influenced by age, gender, race, and body composition. For example, an elderly person could lose half of his/her kidney function before creatinine rises above the upper limit of normal.

Creatinine with eGFR is a more accurate and sensitive assessment of kidney function than creatinine alone. eGFR takes into account the creatinine concentration and other variables including: age, gender, race, and body size.

StatSensor Creatinine calculates eGFR by CKD-EPI, MDRD, and creatinine clearance (CrCl) by Cockcroft-Gault equations.

### eGFR and Chronic Kidney Disease Stages Defined by National Kidney Foundation

Stage	Description	Glomerular filtration (GFR), per 1.73 m <sup>2</sup>
1	Kidney Damage with normal or increased GFR	>90
2	Kidney Damage with mild decrease in GFR	60-89
3a	Moderate decrease in GFR	45-59
3b	Moderate decrease in GFR	30-44
4	Severe decrease in GFR	15-29
5	Kidney failure	<15 (or dialysis)

Kidney Disease: Improving Global Outcomes (KDIGO) CKD Work Group. KDIGO 2012 Clinical Practice Guideline for the Evaluation and Management of Chronic Kidney Disease. *Kidney Int Suppl.* 2013. 3:1-150.

## Another Technology Advance from Nova Biomedical, the World Leader in Whole Blood Testing

Nova Biomedical is the world technology leader in the development of advanced biosensors for whole blood analysis. During the past 36 years, Nova has introduced 20 biosensors, including the industry's first biosensor to directly measure whole blood glucose in 1988 and the industry's first biosensor to directly measure whole blood creatinine in 1996. These biosensors are used routinely in thousands of hospital laboratories and critical care settings around the world in our Stat Profile® series of critical care blood gas analysers. Nova has incorporated its creatinine biosensor technology into a precalibrated, single-use, disposable system that provides rapid screening of kidney function at the point of care.



Nova Whole Blood Measurement Technologies			
Core Technology	Measured Analyte	Abbreviation	Methodology
<b>Amperometric</b> (whole blood or serum)	Creatinine	Crea	Immobilised three enzyme membrane + H <sub>2</sub> O <sub>2</sub>
	Glucose	Glu	Immobilised glucose oxidase membrane + H <sub>2</sub> O <sub>2</sub>
	Glutamate	Glut	Immobilised glutamate oxidase membrane + H <sub>2</sub> O <sub>2</sub>
	Glutamine	Gln	Immobilised glutaminase and glutamate oxidase + H <sub>2</sub> O <sub>2</sub>
	Lactate	Lac	Immobilised lactate oxidase + H <sub>2</sub> O <sub>2</sub>
	Partial Pressure Oxygen	PO <sub>2</sub>	O <sub>2</sub> membrane, O <sub>2</sub> reduction by cathode
<b>Conductivity Potentiometric</b> (whole blood or serum)	Haematocrit	Hct	Electrical resistance, Na <sup>+</sup> corrected
	Acetate	Ace	pH electrode with acetic acid permeable membrane
	Ammonium	NH <sub>4</sub> <sup>+</sup>	Ammonium ISE ionophore in PVC
	Chloride	Cl <sup>-</sup>	Chloride anion ionophore in PVC
	Hydrogen Ion Activity	pH	Hydrogen ion selective glass membrane
	Ionised Calcium	Ca <sup>++</sup>	Calcium ionophore in PVC
	Ionised Magnesium	Mg <sup>++</sup>	Magnesium ionophore in PVC
	Lithium	Li <sup>+</sup>	Lithium ionophore in PVC
	Partial Pressure CO <sub>2</sub>	PCO <sub>2</sub>	pH electrode with CO <sub>2</sub> gas permeable membrane
	Potassium	K <sup>+</sup>	Valinomycin ionophore in PVC
	Sodium	Na <sup>+</sup>	Sodium ion selective glass membrane or sodium ionophore in PVC
	Total Calcium	TCa	Calcium electrode with acidified sample
	Total Carbon Dioxide	TCO <sub>2</sub>	pH electrode, CO <sub>2</sub> membrane, acidified sample
Urea/Urea Nitrogen/BUN	Urea	Urease enzyme membrane and ammonium ISE	
<b>Photometric</b> (lysed whole blood)	Carboxyhaemoglobin	COHb	Multi-wavelength spectral scanning of haemolyzed RBCs
	Deoxyhaemoglobin	HHb	Multi-wavelength spectral scanning of haemolyzed RBCs
	Methaemoglobin	MetHb	Multi-wavelength spectral scanning of haemolyzed RBCs
	Oxygen Saturation	SO <sub>2</sub> %	Multi-wavelength spectral scanning of haemolyzed RBCs
	Oxyhaemoglobin	O <sub>2</sub> Hb	Multi-wavelength spectral scanning of haemolyzed RBCs
	Sulphaemoglobin	sHb	Multi-wavelength spectral scanning of haemolyzed RBCs
	Total Haemoglobin	tHb	Multi-wavelength spectral scanning of tt RBCs
<b>Photometric</b> (whole blood, nondestructive)	Cell Density	CD	Multi-wavelength fiber optic reflectance
	Haemoglobin	Hb	Multi-wavelength fiber optic reflectance plus conductivity, sodium correction
	Oxygen Saturation	SO <sub>2</sub> %	Multi-wavelength fiber optic reflectance (oximetry)

# Simple, Fast, Accurate Creatinine/eGFR Testing at the Point of Care



## Rapid, Easy-to-Use, Point-of-Care Creatinine Testing

### No User Calibration

StatSensor biosensors are precalibrated and immediately ready to use.

- No analyser or biosensor preanalytical steps, such as using a calibrating reagent or entering a calibration code in the analyser, are necessary.
- By eliminating this calibration step, StatSensor is actually simpler to operate than most meters used by diabetics at home to self-test for glucose.

### Simple Test Procedure

A bright, easy-to-read, colour screen prompts the user through the StatSensor test procedure.

1. Place the sensor in meter
2. Stick finger to create blood drop
3. Apply sensor to blood drop
4. Read creatinine results



### Creatinine Results in 30 Seconds

Assessment of kidney function can be made in real time at the point of care. Patient care treatments that require prior kidney assessment can be expedited.

### Creatinine with eGFR Enhances Renal Function Assessment

StatSensor can calculate eGFR by three different equations. eGFR and creatinine results are displayed along with normal ranges. Abnormal results are prominently flagged with colour highlighting, text, and symbols.

### Cleared for Use by Point-of-Care Personnel

StatSensor can be used by nurses, radiology technicians, or other non-laboratory personnel.

### Capillary Blood Sampling is Virtually Painless

StatSensor uses a very small, 1.2 microlitre blood sample (less than one drop). Samples can be obtained with virtually no pain from a capillary finger stick using the same type of lancet used by diabetics at home to self-test for glucose. No venous blood drawing or phlebotomist is needed.



## StatSensor Creatinine is Comparable to Hospital Laboratory Testing

StatSensor Creatinine's whole blood biosensors provide excellent correlation to central laboratory plasma reference methods. StatSensor accuracy combined with fast analysis in 30 seconds and calculated eGFR enable rapid screening of kidney function at the point of care.

## StatSensor shown to be more accurate than the laboratory method in a large medical center study

In a 670 patient study funded by the International Society of Nephrology, the South Africa Medical Research Council and the University of Witwatersrand, Johannesburg, South Africa, the Nova POC StatSensor Creatinine/eGFR meter was more accurate than the central laboratory IDMS-traceable Jaffe methodology in estimating GFR when both methods were compared to MEASURED GFR (iohexol).<sup>1</sup>


## StatSensor Measures and Eliminates Interference Due to Haematocrit

Measurements made on whole blood samples are subject to interference due to varying haematocrit levels. StatSensor measures and corrects for varying haematocrit levels. Accurate results are obtained throughout a broad range of haematocrit.


## Interfaces with Hospital and Laboratory Information Systems

NovaNet provides comprehensive connectivity, management, and control for StatSensor POC testing. This system provides the critical components for successful POC testing including:

- LIS/HIS interface and connectivity management
- Management of patient and QC data
- Regulatory compliance via audit logs, reports, and records

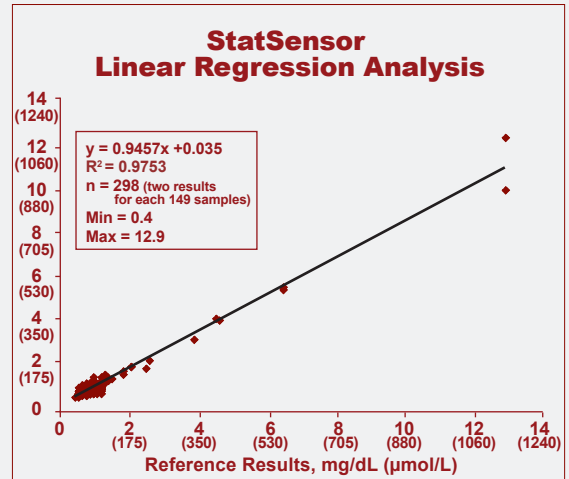


StatSensor Whole Blood vs Central Laboratory Plasma Reference Method (Enzymatic)

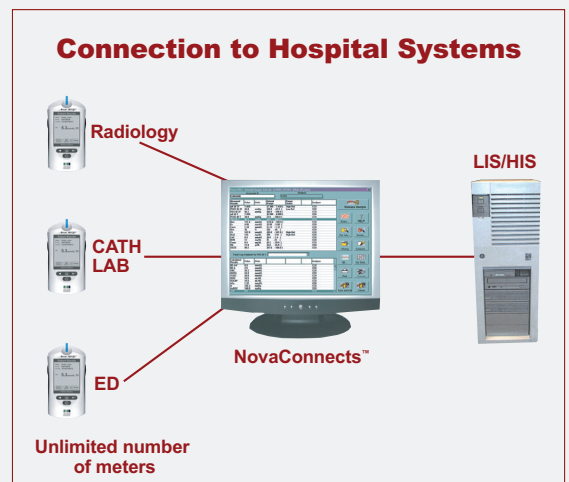

1. George J et al. Evaluating chronic kidney disease in rural South Africa: comparing estimated glomerular filtration rate using point of care to iohexol measured GFR. CCLM 2021

# Point-of-Care Creatinine Testing for Radiology and Cardiac Labs

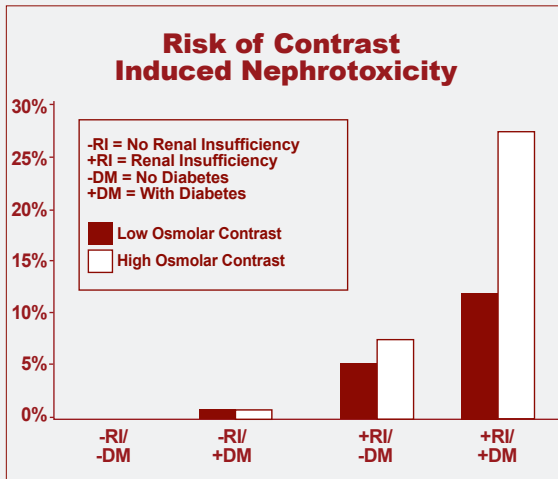

The primary risk factor for contrast induced nephrotoxicity is pre-existing renal dysfunction, especially diabetic nephropathy.<sup>2</sup>

## For Radiology and Cardiac Catheterization

### Minimise Risk of Contrast Induced Nephropathy

Contrast induced nephropathy (CIN) is the third most common cause of acute renal failure in hospitalised patients.<sup>1</sup> CIN is associated with prolonged hospital stay, adverse cardiac events, and high mortality. An increasing number of imaging procedures require the use of intravenous contrast, and the patient population subjected to these procedures is progressively older and has more pre-existing conditions. The benefits of a fast, easy POC creatinine assay are:

- StatSensor can minimise the risk of CIN by providing comprehensive assessment of renal impairment prior to contrast imaging.
- Comprehensive renal assessment allows renal protective strategies to be considered.
- Testing can be easily performed by radiology personnel.

### Improve Productivity and Workflow

Radiology departments and cath labs are highly specialised and costly to operate. Procedures are carefully scheduled in order to maintain productivity of personnel and procedure rooms. Patients arriving for contrast imaging procedures without valid creatinine results must be sent to the laboratory for creatinine testing. Obtaining creatinine/eGFR results from the lab can delay the procedure for hours or require rescheduling for another day. Productivity is lost when schedule openings occur.

- StatSensor POC creatinine/eGFR testing can prevent costly procedure room openings, and loss of personnel and equipment productivity.
- A simple, inexpensive StatSensor test can prevent cancellation or rescheduling of a scan when kidney function testing is needed.

### Improve Patient Satisfaction

Long delays or even cancellation and rescheduling of procedures can occur when renal status is unknown at the patient's time of appointment. Patients may be forced to wait for hours for their creatinine/eGFR blood test or to come back another day. Patient dissatisfaction is likely if these events occur.

- StatSensor Creatinine provides rapid assessment of renal function in 30 seconds with virtually no pain from a finger stick.
- StatSensor POC testing in the department can prevent patient dissatisfaction when renal function must be determined.

1. Hou, S. H., Bushinsky, D. A., Wish, J. B., Cohen, J. J., & Harrington, J. T. (1983). Hospital-acquired renal insufficiency: A prospective study. *Am J Med*, 74(2), 243-248.

2. Rudnick, M. R., Goldfarb, S., Wexler, L., Ludbrook, P. A., Murphy, M. J., Halpern, E. F., ... VanFossen, D. B. (1995). Nephrotoxicity of ionic and nonionic contrast media in 1196 patients: A randomised trial. *The Iohexol Cooperative Study. Kidney Int*, 47(1), 254-261.

# Point-of-Care Creatinine Testing for Oncology, Emergency Department, and ICU

## For Oncology

### Reduce Medication Risks

Renal impairment is frequent in cancer patients, and these patients are at high risk of drug-induced renal toxicity. This has implications for drug selection and dosing. Many chemotherapy drugs are prescribed close to the maximum therapeutic dose. Renal function determines whether there is a need to use the drug more sparingly or to avoid the drug entirely. This is especially true for drugs that are cleared primarily by the kidney and for drugs with established nephrotoxicity.

- StatSensor provides a rapid, accurate assessment in 30 seconds of renal function including eGFR.
- StatSensor testing can be easily performed in the oncology clinic.
- Renal safety for chemotherapy patients is ensured.

### Improve Patient Satisfaction

Many chemotherapy outpatients need to have their kidney function assessed before their chemotherapy session. Obtaining a venous blood sample and sending it to a laboratory for creatinine/eGFR testing could delay chemotherapy for several hours. A chemotherapy session that might take a couple of hours could turn into a full day.

- With StatSensor Creatinine, a capillary blood sample can be obtained and renal function can be quickly determined.
- Treatment delays and patient dissatisfaction are avoided.

## For the Emergency Department

### Improve Triage, Expedite Patient Care


Various patient care protocols such as chest pain and stroke protocols may require rapid assessment of renal function in order to triage the patient for radiologic procedures. Many other ED protocols require timely administration of drugs and adjustment of drug dosage based on renal function.

StatSensor Creatinine accurately assesses renal function with a simple test that can be easily performed in the ED in 30 seconds.

## For the ICU

Acute renal failure (ARF) is a common occurrence among patients admitted to the ICU. Results of a multinational study suggest that ARF requiring renal replacement therapy occurs in 5% to 6% of ICU patients.<sup>4</sup> Mortality for ARF in the ICU setting may be as high as 70%. Early detection and therapy is critical in reducing morbidity and mortality from ARF.

- StatSensor Creatinine with eGFR provides real time assessment of kidney function with a simple test in 30 seconds.
- Testing can be performed at the bedside by nursing or other POC test personnel.

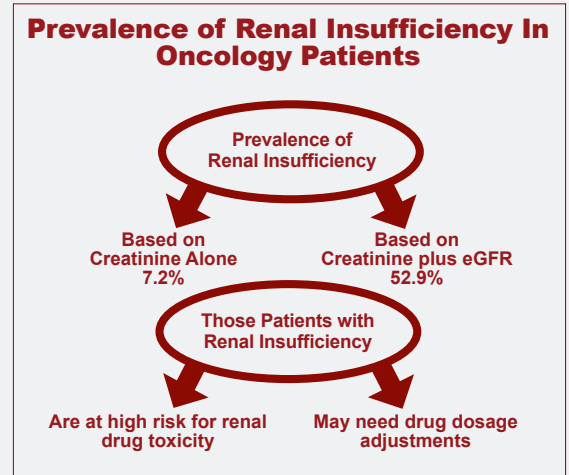

*The Insuffisance Renale et Medicaments Anticancereux (IRMA) study of more than 4,600 cancer patients found that renal insufficiency is common among cancer patients—a problem that complicates drug selection and dosing.<sup>3</sup>*

Patient Result		Op: 8636
Name: A Unit-1		⏪ 13:42
ID: 194546646	Strip Lot: 0206806099	
Darryl Smith		
<b>Creat</b>	<b>97</b> $\mu\text{mol/L}$	Normal 53-106
<b>GFR</b>	<b>54</b> mL/min/1.73 m <sup>2</sup>	
⏪ Reject	OK	Comment ⏩

3. Launay-Vacher, V., Oudard, S., Janus, N., Le Tourneau, C., Rixe, O., Pourrant, X., ... Morere, J. (2006). Renal insufficiency in cancer patients: Prevalence and implications on anticancer drugs management. Results of the IRMA study [Abstract]. Retrieved from [http://meeting.ascopubs.org/cgi/content/abstract/24/18\\_suppl/8603?sid=8b13b8a8-8eab-4e9e-a46e-f5437fb1356b](http://meeting.ascopubs.org/cgi/content/abstract/24/18_suppl/8603?sid=8b13b8a8-8eab-4e9e-a46e-f5437fb1356b)

4. Uchino, S., Kellum, J. A., Bellomo, R., Doig, G. S., Morimatsu, H., Morgera, S., ... Ronco, C. (2005). Acute renal failure in critically ill patients: A multinational, multicenter study. *JAMA*, 297, 813-818. doi:10.1001/jama.294.7.813

# StatSensor® CREAT



## StatSensor Specifications

Test Measured: ..... Creatinine  
 Tests Reported: ..... Creatinine, eGFR  
 Test Time: ..... 30 seconds  
 Test Strip Volume: ..... 1.2 µL  
 Test Methodology: ..... Electrochemistry  
 Weight: ..... 266g (0.6 lbs)  
 Size: ..... 153 mm x 82.5 mm x 46 mm  
 (6.0 in x 3.25 in x 1.8 in)

### Sample Types and Operating Modes:

Whole Blood: ..... Arterial, Venous, Capillary

### Measurement Range:

Creatinine: ..... 27 – 1056 µmol/L (0.3 – 12.0 mg/dL)

### Operating Ranges:

Temperature: ..... 15°C – 40°C (59°F – 104°F)  
 Altitude: ..... Up to 4,500 meters (15,000 feet)  
 Humidity: ..... 10% to 90% relative humidity

### eGFR and CrCl Equations:

CKD-EPI, MDRD, Cockcroft-Gault

### Certifications and Compliance:

Nova Biomedical complies with FDA Quality System Regulations and is certified to EN ISO 13485:2016. StatSensor meters are tested according to EN 61010-1:2010 EN 61010-2-101:2015

### FDA Labeling:

For in-vitro diagnostic use

### Data Storage:

Patient Tests: ..... 1,000  
 QC Tests: ..... 200  
 Operators: ..... 8,000

### Connectivity:

Meter Data Output: ..... RJ-45 Ethernet Port  
 Protocol: ..... TCP/IP Ethernet 100 Mbit/1 Gbit  
 Standard: ..... POCT1-A Compliant  
 Setup Program: ..... NovaNet™ Web-Based  
 Instrument Manager Software

### Battery Information:

Type: ..... 3.7V Li Polymer Rechargeable Battery

## StatSensor Xpress Specifications

Tests Measured Blood: ..... Creatinine  
 Creatinine Methodology: ..... Enzyme, Amperometric Creatinine  
 Test Results: ..... mg/dL, µmol/L  
 Sample Type Whole Blood: ..... Capillary, Arterial, and Venous  
 Creatinine Test Range: ..... 0.30 – 12.0 mg/dL  
 27 – 1056 µmol/L  
 Test Time: ..... 30 seconds  
 Test Strip Volumes: ..... 1.2 µL  
 Memory storage: ..... 400 Tests  
 Battery Life (nominal): ..... 600 Tests  
 Battery Type: ..... 3V coin cell  
 Data Cable Serial or USB

## StatSensor Test Strips

Measurement Range: ..... 27 – 1056 µmol/L  
 0.3 – 12.0 mg/dL  
 Test Strip Stability: ..... Room temperature 3 months  
 Refrigerated 12 months  
 Strip Packaging: ..... Box of 50 test strips  
 2 vials of 25 test strips



Nova Biomedical Patent Numbers: US 6,767,441 CA 2,470,465

Nova Biomedical® and StatSensor® are a registered trademark of Nova Biomedical Corporation.

**nova**  
biomedical  
novabiomedical.com



**Nova Biomedical Headquarters:** 200 Prospect St., Waltham, MA 02454 U.S.A., +1-781-894-0800 800-458-5813, FAX: +1-781-894-5915, Int'l FAX: +1-781-899-0417, e-mail: info@novabio.com  
**Nova Biomedical Australia ANZ Pty. Ltd.:** 5/372 Eastern Valley Way, Chatswood, NSW, 2067, Australia, TEL: +61(0) 2 9417 0193, email: AU-info@novabio.com  
**Nova Biomedical Benelux B.V.:** Korenmoelen 22, 5281 PB, Boxtel, The Netherlands, TEL: +31(0)733032701, e-mail: benelux-info@novabio.com  
**Nova Biomedical Brasil:** Rua Massena, 107, Jardim Canadá, Nova Lima - MG, CEP: 34007-746 Brasil, TEL: +55-31-3360-2500, e-mail: BR-info@novabio.com  
**Nova Biomedical Canada, Ltd:** 17 – 2900 Argenta Road, Mississauga, Ontario L5N 7X9 Canada, TEL: +1-905-567-7700 800-263-5999, FAX: +1-905-567-5496, e-mail: CA-info@novabio.com  
**Nova Biomedical France:** Parc Technopolis - Bât. Sigma 3 Avenue du Canada 91940 Les Ulis Courtaboeuf, France, TEL: +33-1-64 86 11 74, FAX: +33-1-64 46 24 03, e-mail: FR-info@novabio.com  
**Nova Biomedical GmbH, Deutschland:** Hessenring 13 A, Geb. G, 64546 Mörfelden-Walldorf, Germany, TEL: +49-6105 4505-0, FAX: +49-6105 4505-37, e-mail: DE-info@novabio.com  
**Nova Biomedical Iberia, S.L.:** c/Vic 17, Planta 3A 08173 Sant Cugat del Vallès, Barcelona, Spain, TEL: +34 935531173, e-mail: ES-info@novabio.com or PT-info@novabio.com  
**Nova Biomedical Italia S.r.l.:** via Como, 19 - 20045 Lainate (MI), Italy, Tel: +39 02 87070041, Fax: +39 02 87071482, e-mail: IT-info@novabio.com  
**Nova Biomedical K.K., Japan:** Harumi Island Triton Square Office Tower X 7F, 1-8-10 Harumi, Chuo-ku, Tokyo 104-6007, Japan, TEL: 03-5144-4144, FAX: 03-5144-4177, e-mail: jp-info@novabio.com  
**Nova Biomedical New Zealand:** Regus Constellation Drive - Candida Building 4 Level 3/61 Constellation Drive, Rosedale, Auckland, 0630, New Zealand, TEL: +61(0) 2 9417 0193, email: AU-info@novabio.com  
**Nova Biomedical Schweiz GmbH:** Herostrasse 7, 8048 Zürich, Switzerland, TEL: +41-41-521-6655, FAX: +41-41-521-6656, e-mail: CH-info@novabio.com  
**Nova Biomedical U.K.:** Innovation House, Aston Lane South, Runcorn, Cheshire WA7 3FY United Kingdom, TEL: +44-1928 704040, FAX: +44-1928 796792, e-mail: UK-info@novabio.com

Specifications current as of revision date.

177E V4 INT EN 10/10/24