

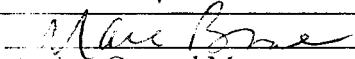


ADDENDUM TO CERTIFICATE OF COMPLIANCE Page 2/2
RELEASE TEST RESULTS

Customer:	NIPRO DIAGNOSTICS INC.
Description:	TEST STRIPS, NIPRO STATSTRIP™ GLUCOSE
Catalog Number:	47940
Nova Biomedical P/N	47940
Lot No.	2125310309

(Spiked whole blood) Target Range:30-50		
Mean (mg/dL) N=100	SD (mg/dL) ≤4	YSI (reference, mg/dL)
45	2.9	47.5
(Spiked whole blood) Target Range:51-110		
Mean (mg/dL) N=100	CV% ≤6.0	YSI (reference, mg/dL)
96	4.7	95.6
(Spiked whole blood) Target Range:111-150		
Mean (mg/dL) N=100	CV% ≤5.0	YSI (reference, mg/dL)
130	3.4	132.8
(Spiked whole blood) Target Range:151-250		
Mean (mg/dL) N=100	CV% ≤5.0	YSI (reference, mg/dL)
218	3.2	218.0
(Spiked whole blood) Target Range:251-400		
Mean (mg/dL) N=100	CV% ≤5.0	YSI (reference, mg/dL)
357	3.7	346.3

Accuracy Criteria 1: Total within +/-12 mg/dL or +/- 12% of reference method ≥ 95%	Accuracy Criteria 2: Total within +/-15 mg/dL or +/- 15% of reference method ≥ 98%
100%	100%


 Quality Control Manager
 Nova Biomedical Corporation
 Date: 11/11/25
 200 Prospect Street Waltham, MA 02454